

Provetscan SR-2C

High resolution ultrasound scanner
for cattle and horses



Manual User

Introduction

The SR-2C wireless probes are cutting-edge ultrasound instruments for reproductive diagnosis.

Unlike traditional veterinarian ultrasounds with a cable connected to the probe, the SR-2C is wireless.

The SR-2C probe has a built-in ultrasound image processor, power source and wireless signal module to connect to the main unit.

Now, the main unit is replaced by a tablet. The probe acts as a WiFi 2,4 and 5 GHz access point and shows the image through the VetUS application.

The purpose of this manual is to provide a general overview. Read it carefully before using the device.

1. Features

Screen: iPad series

Scan: Linear

Probe: 6.5 / 8 MHz transducer

Depth: 40-100 mm

Display Mode: B Mode, BM, Colour, PDI and PW

Grey scale: 256 levels

Battery life: 6 hours of continuous work, or 14 hours on stand-by.

Size of main unit: 180mm x 63mm

Examination wire length: 1.8 m.

Weight: 400 g

2. Start

For your protection, please read these safety instructions before starting up the probe.

2.1 Unpackaging

The SR-2C probe is carefully packaged to avoid damage during transport. Before opening the case, look for any visible damage on the outside that may have occurred during transport.

Items must be verified to ensure that they were all received. The table below shows the items that must be included:

Items	Included
SR-2C Ultrasound probe	✓
USB cable for charging	✓
Transport suitcase	✓
Wireless charger	✓

Each item must be examined to report defects or damage that may have occurred during shipping, despite careful packaging. If this occurs, please immediately contact your distributor to report the issue.

2.2 Installation the app

If the VetUS application is not installed on your tablet, please go to the App Store and search for it. Installation is free of cost.

2.3 Connecting the probe

The WiFi connection and battery life indicators must be turned off before connecting the probe.

Press the **(1)** ON/OFF button. The battery life indicator **(5)** tells you how much batter you have left. This indicator is split into four levels.

Seconds after connecting the probe, the **(6)** WiFi connection indicator will blink, showing it is ready to be connected to the mobile device.

Press the **(1)** ON/OFF button for three seconds to turn the probe off. When the probe is turned off, the light-up indicators **(5)** and **(6)** will also be turned off.

Additionally, this WiFi probe has another 3 buttons: **(2)** to modify depth, **(3)** to lower gain and **(4)** to raise gain.



2.4 WiFi Connection

When the probe is waiting for a WiFi connection, go to iPad settings and search for the probe's SSID. The SSID is, for example, "SR-1C GVCJDA001." Connect the SSID with the password, which is the serial number engraved on the bottom of the probe. The password would be, for example, "WRPCJDA001."


The password's letters must be input as small letters, not capital.

Once connected to the WiFi network, start up the VetUS application. When we begin working with the probe, the WiFi signal indicator will no longer blink the same way.

3. App functions


Once the probe is connected, start up the app and the main screen will appear.
No image is shown when you open the app for the first time.

3.1 Ultrasound scan

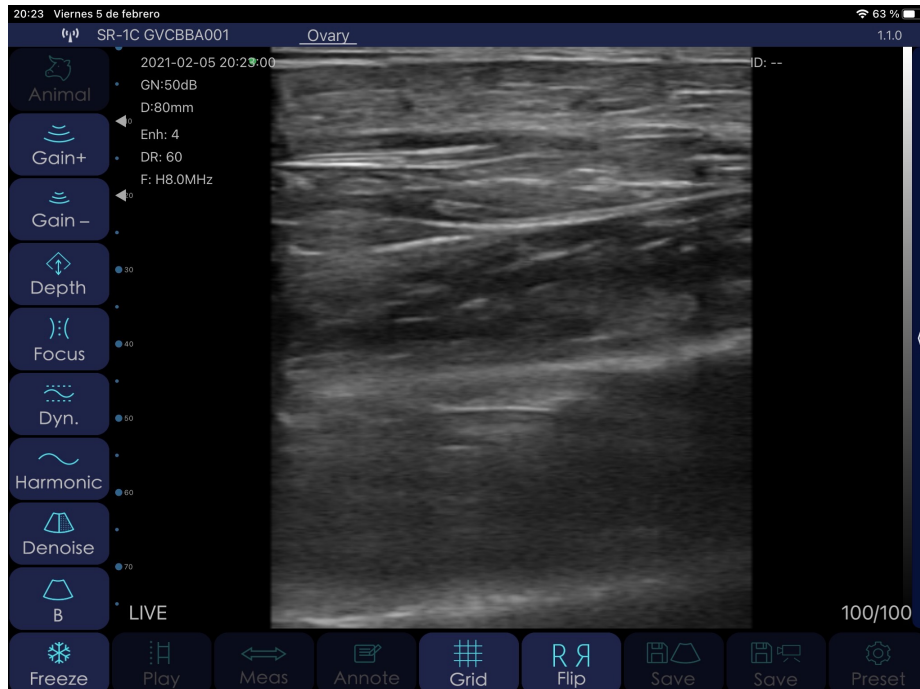


The connected probe's network is shown on the upper left.

Description of app buttons:

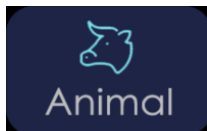

When you press the "Animal" button, a window pops up to enter the type of animal and the following fields:

When you press on the Animal image, three possibilities are displayed: bovine livestock, equine livestock and unspecified.

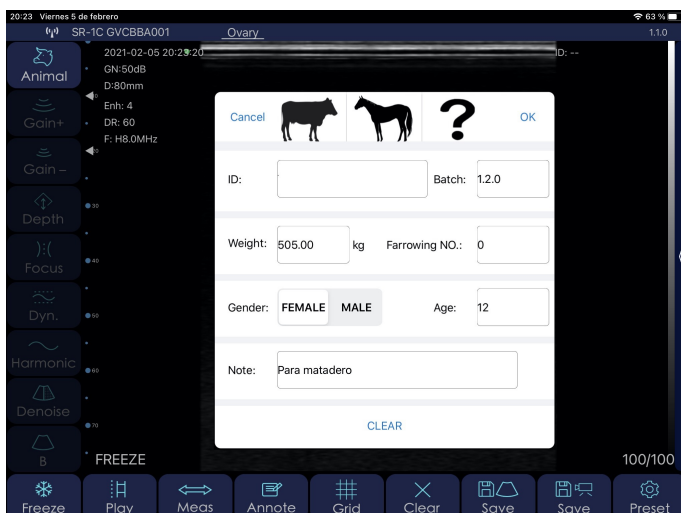

ID: ID number of animal with alphanumeric characters.

Batch: Batch number with alphanumeric characters.

Weight: Animal weight.

Farrowing No.: farrowing num. (not shown when Male is selected)

Gender: Select male or female

Age: Animal age

Note: To insert notes or diagnoses on the animal.

All this fields are recorded in the CSV file generated when you export the data stored in the application.

GAIN is the ability to modify the ultrasound's amplitude (magnitude of the ultrasound), changing the brightness of the image. General gain changes affect the entire image equally.




Raises scan gain up to 105 dB.



Lowers scan gain down to 30 dB.

Depth


You can work with greater or lesser depth. Press the following button.

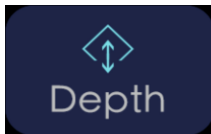

Selects the 4 zooms available: 40, 60, 80 and 100 mm.
You can also adjust this by swiping your finger in the middle of the screen, up or down.

Focus


This is the maximum concentration of lines inside the ultrasound beam. This means you must place the focus in the area of interest.

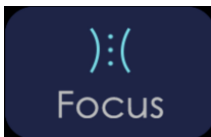

Adjusts the focus to obtain the best image, raising or lowering it on the screen.

Dinamic range

This process raises or reduces the difference between ultrasounds of greater or lesser amplitude. With the Dynamic Range, you vary the amplitude of the signal shown in grey scale. At a low setting, the image has greater contrast. When increased, the image is softened.

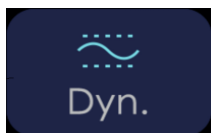

The value range goes from 40 to 110. Lower settings increase the image's contrast.

Frequency

You must select the suitable frequency for the examination. With greater frequency, there is lesser penetration, and vice-versa. With the Harmonic Image, the return signal not only includes the base frequency, but also the frequency matching its double. The ultrasound processes both signals separately, eliminating the basic frequency with its noise.



Changes the scan frequency from 6.5 to 8 MHz

Denoise


It's an algorithm to eliminate spotting, to soften the image of the organs and to make the edge of tissue softer. Lower settings increase image noise, and high settings reduce it.



Reducing image grain. Range of value from 0 to 3.

B Mode

Two-dimensional representation of the ultrasounds, shown as light-up dots with variable clarity, depending on frequency and depth.

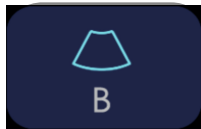

Press this button to show the ultrasound options with modes BM, Colour, PDI and PW. More information in section 3.2.

Freeze

Put the probe on stand-by by freezing the image to examine aspects more thoroughly. The probe's consumption is less than half of LIVE mode consumption.



Freeze and unfreeze (live) the image.

Play


This equipment permanently records a video sequence (Cine Loop) of the last 100, 200, 500 or 1000 frames. As new photograms are recorded, the previous ones are deleted. To change the Cine Loop's size, go to the settings button.



Once the image is frozen, this plays the saved Cine Loop.

Measurements


You have three measurement modes: length, area/circumference, and contour line. Once the method is selected, take the measurement with your finger on the screen, placing the cursors where you wish to measure for length mode, or by drawing the desired area for area, circumference, and contour line modes.

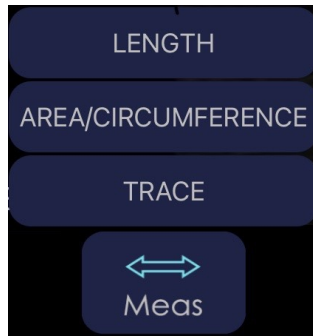


Press this button to drop-down the menu that shows length, area/ circumference and contour line measurements, and measurements to calculate gestation if you selected bovine or equine livestock.

Gestation age calculation:

If you have selected bovine livestock, you make take the following measurements:

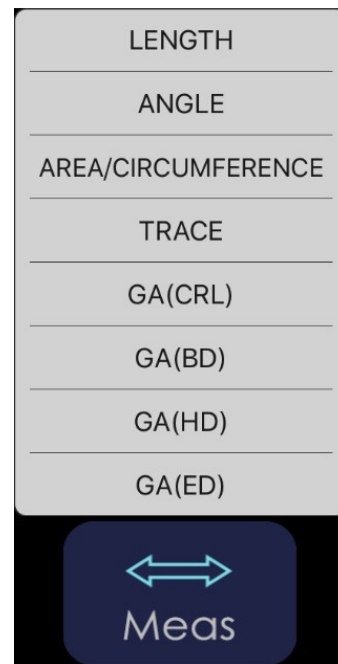

GA (CRL): Crown to Rump Length. This calculates the gestation age based on the distance from the crown of the head to the rump of the foetus.

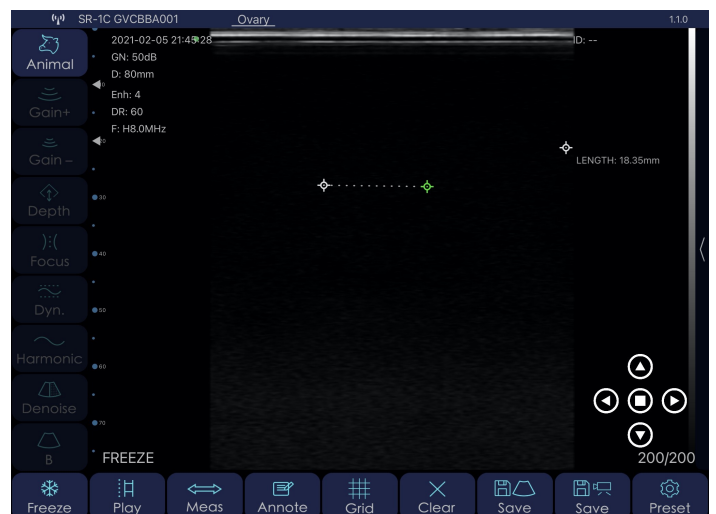

GA (BD): Body Diameter. Calculates the gestation age based on the maximum diameter of the foetus' body/ thorax.




GA (HD): Head Diameter. Calculates the gestation age based on the maximum diameter of the foetus' head.



GA (ED): Eye Diameter. Calculates the gestation age based on the maximum diameter of the foetus' eye.

If you have selected equine livestock, you make take the following measurements:


Note: once the gestation measurement has been taken down, the estimated delivery date (EDD) is automatically shown

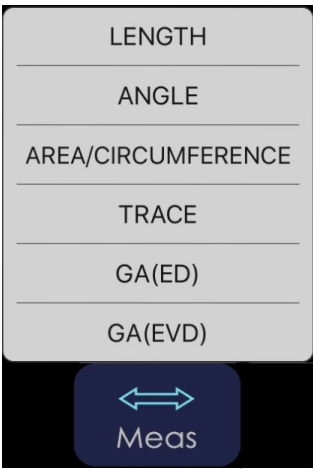

Meas Ar



GA (ED): Eye Diameter. Calculates the gestation age based on the maximum diameter of the foetus' eye.

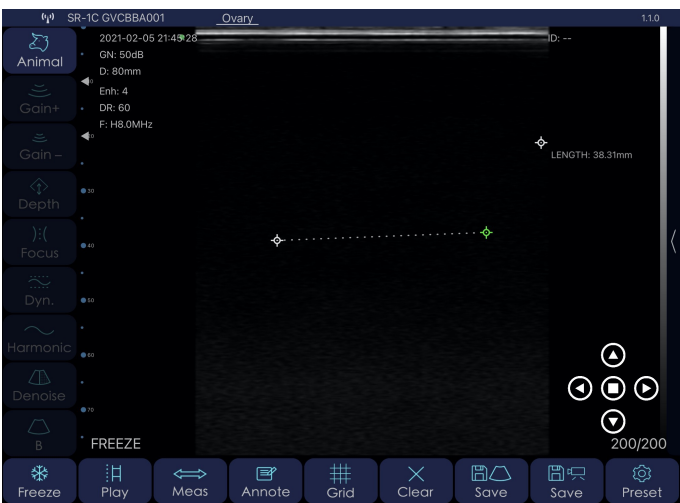

GA (EVD): Embryonic Vesicle Diameter. Calculates the gestation age based on the maximum diameter of the embryonic vesicle.

Annotations

You can make notes anywhere on the screen. You may input alphanumeric characters. If, after leaving a comment on the screen, you save the image, then the note will be included.



This creates a window at the location on the screen where you pressed with your finger, and the keyboard pops up to write.


Grid To have orientation of the measurements organs may have, the size of the foetus, or other structures, you can work with a grid superimposed on the screen.



Activates or deactivates the grid on the screen.

Change scan direction

In scan mode, press this button to change the scan direction.

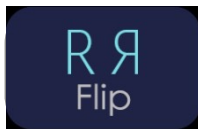

A green dot on the upper right or left of the ultrasound image shows this. If the dot is to the right, it means the image will be viewed starting from the right as you insert the probe into the rectum. If the dot is on the left, it is the other way around.

Clear



Press this button to clear all measurements and annotations on the screen.

Save Image and diagnostic data

Once the image is frozen, you have saved the 100, 200, 500 or 1000 images in the Cine Loop that may be saved in the film roll. Swipe your fingers in the middle of the screen to the left or the right to select the image you wish to save.




Save a screen shot of the image selected from the film roll.

At the same time an image is saved, a window pops up showing the saved diagnostic data. This window automatically closes after about three seconds.




If you wish to see the table with all diagnostic data, press the button the upper-right part of the image:


An example of diagnostic data saved in the application's memory might be:

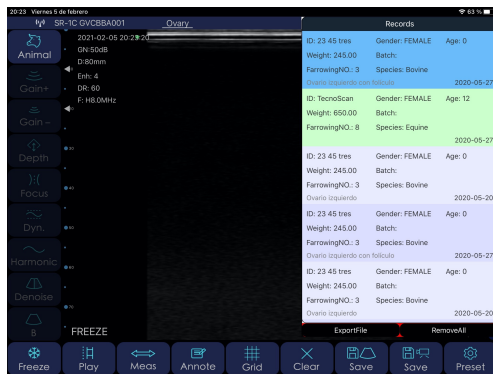


This archive will be included inside a folder with today's date.

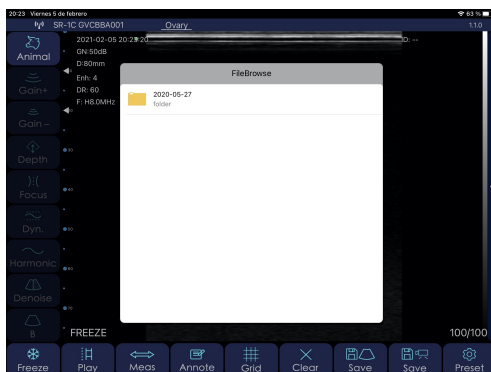


You can preview the CSV file by pressing on the file icon.

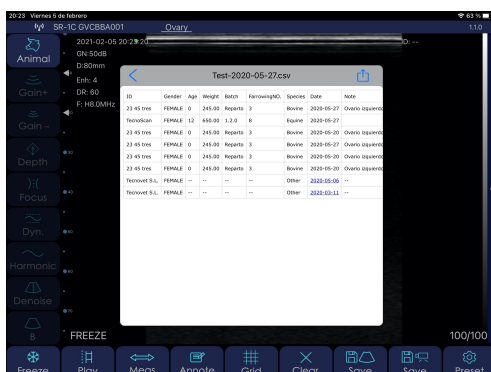


Once this table is displayed, you have the possibility of deleting all diagnostic data (RemoveAll) or exporting these data (ExportFile). If you select "ExportFile," a new window pops up to name the CSV file generated, along with today's date:

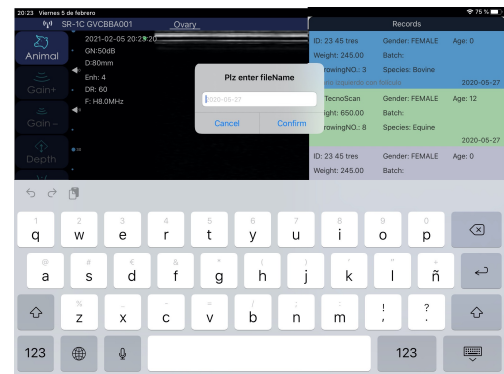

It is here inside this folder where you can find the CSV file. This file automatically opens with Excel, although you can also open it in Numbers.

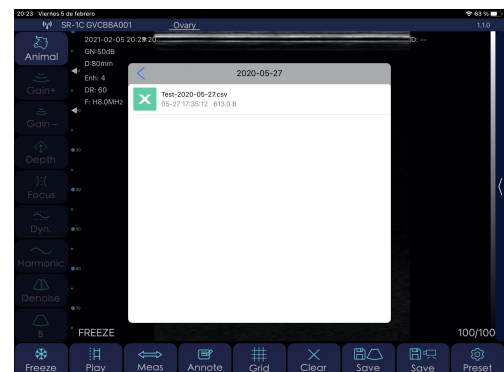

You can now share this file with applications on your iPad and save them in these applications with the Files application:

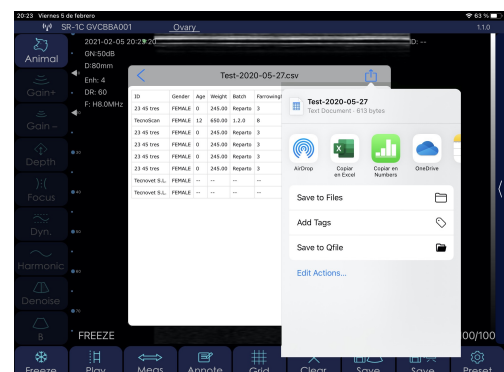

If you wish to access the saved folders holding the CSV files, press on the following icon, on the upper right-hand part of the screen:

Save Video

When you freeze the image, like in the previous case, the Cine Loop stops recording, and the entire photogram sequence can be saved. Videos will last 10, 20, 50 or 100 seconds.



This saves in the Cine Loop you were storing until you pressed Freeze in the film roll.


Settings

Shows the settings you can adjust on the equipment.




Press this button for a pop-up window to modify the following parameters:

- Show information: shows or hides information on screen parameters.
- RightHand: If you activate this, the buttons that are shown by default on the left will be shown on the right. This is more comfortable for users who conduct the examination with their left hand.
- Compound: Processes the image more so it appears softer, but the image is slower.
- Cine Loop: temporary recording made as you conduct the ultrasound. The last 100, 200, 500 or 1000 photograms will be saved.
- Auto Freeze: Selection of the waiting time from which the image freezes to save energy.
- Wireless Channel: Adjusting WiFi, 2.4 GHz or 5 GHz channels. In the event of conflict with other WiFi network channels nearby, the channel can be changed to prevent interference. If you change the channel, you will lose the probe password and you will need to enter it again. See point 2.4 of this manual.

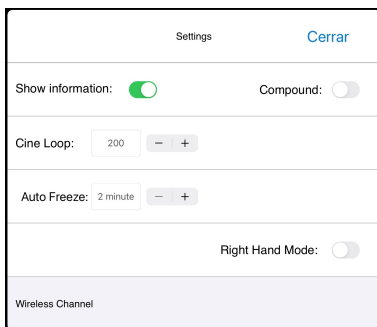

In addition to these buttons, there is a drop-down menu on the upper left where you can select up to 6 pre-settings.

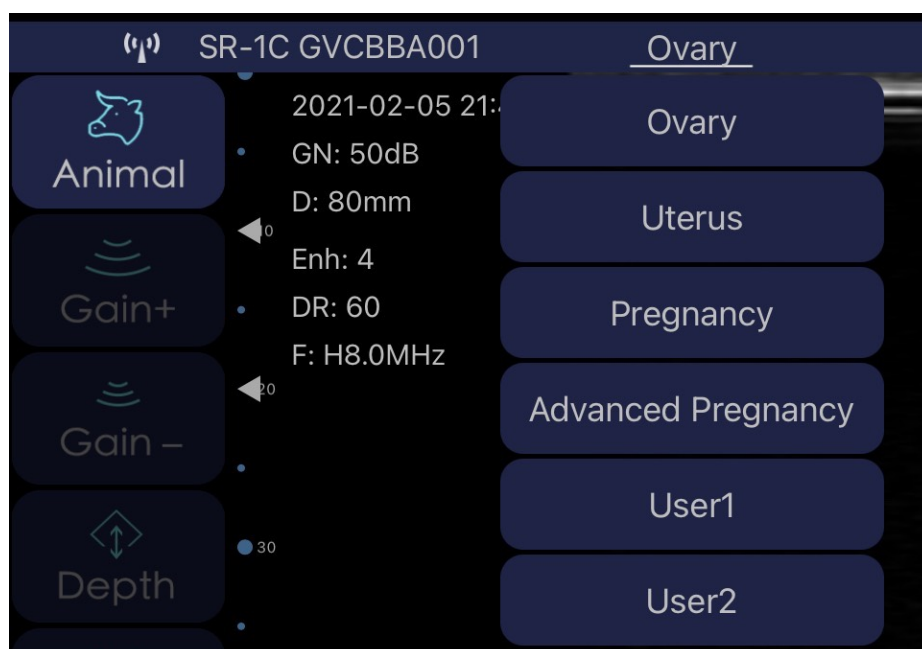

The first 4, **OVARY**, **UTERUS**, **PREGNANCY** and **ADVANCED PREGNANCY** are factory settings. Although you can modify them, every time you press one of them, you will return to the memorised values.

The next 2, **USER1** and **USER2**, can be adjusted by users. The last values right before closing the app are always memorised, so that the latest settings are shown when you return to work.




Finally, if you press the vertical bar on the right-hand side, you can drop down the 8 sections of the **TGC** to adjust gain as desired for each zone.

If you slide each one of the buttons to the right, you increase gain at that screen height. Sliding it to the left decreases gain. If you press **RESET**, you return the factory TGC settings.

Note: Only buttons lit up in white are available. On the lower right, you can see a number that shows the current image recorded in the Cine Loop (100/100). The user can swipe their finger in the middle of the screen to the left or the right to select the desired image. This means that any photogram, from 1/100 to 100/100, may be saved by pressing the button "Save image."


The lower left of the screen will read "LIVE" when the probe is in motion or will show "FREEZE" when the image is frozen.

3.2 Special functions: Settings in BM and Doppler Mode

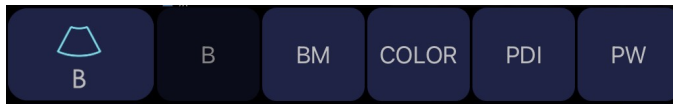

As explained in point 3.1, if you press the “B” button, you access the rest of the ultrasound machine’s special functions: BM, Colour, PDI and PW modes.


3.2.1 BM Mode

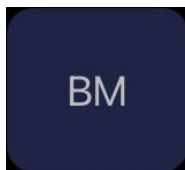

COLOR


The M Mode shows the movement of all information collected along a line we place over the image. This provides for a graphic analysis of surfaces that are in movement, and to quantify them. This button combines the B (2D ultrasound) and M (in movement) modes. M signals are shown on the monitor’s horizontal axis.

3.2.2 COLOUR Mode

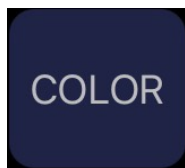

The Doppler effect is the change in frequency felt by a given transferred wave that is detected by a receptor, due to certain conditions in movement between the receptor, the transmitter, the medium and the reflective surfaces. The change in frequency is caused by the effect producing movement on the speed of the wave between the transmitter and the receptor.

The colour Doppler encodes the flow direction in two colours. Conventionally, red is the flow approaching the transducer, while blue is the one moving away from the transducer. Changes in tone also indicate the circulatory speed: bright tones show high speed, while dull tones show a low circulatory speed. When flow speed is too high, the laminar flow grows turbulent (multiple speeds in one same location at one same time). In this case, it is represented by adding green to the underlying red-yellow or blue-white, thus creating a mosaic pattern.




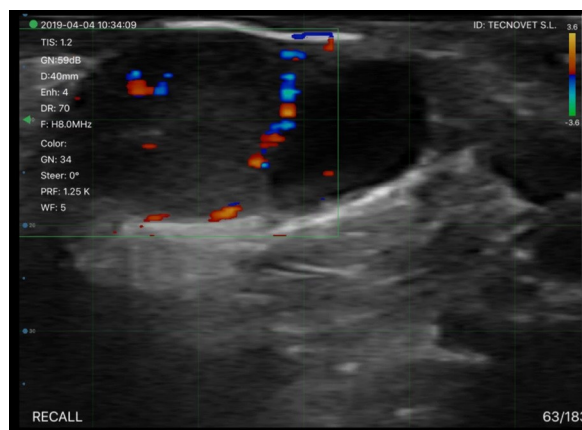


Colour mode settings:



Steering consists of tilting the Doppler image. The angle of incidence between the propagation line and the flow direction must be 45° - 60° . Given the formula that governs the Doppler principle, small changes in angles greater than 60° lead to speed errors which are much greater than small changes to angles less than 60° .

In this case, the transducer must be relocated, or the Steering technique used. You can vary the angle at $\pm 7^{\circ}$.

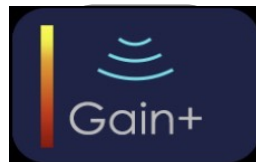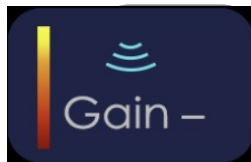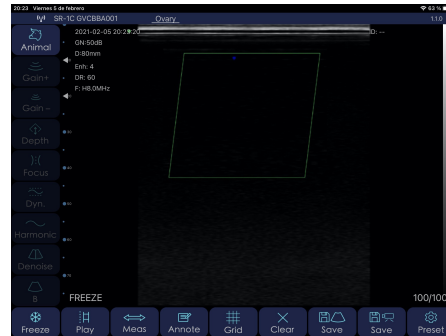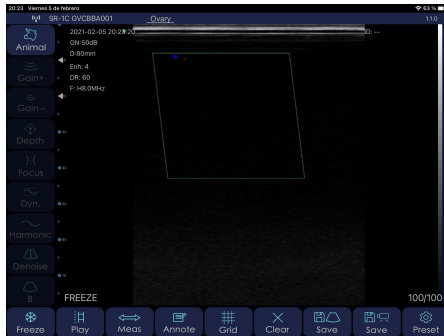


In addition to adjusting the B-mode gain, you may adjust Gain in Colour mode independently. This amplifies or decreases the signal received, seeing more or fewer blood vessels. It is advisable that the colour fill the vessel with light, from one side to the other.




The Pulses or Pulse Repetition Frequency is a parameter set by the operator. It defines the rate or frequency with which the transducer will emit ultrasound pulses to show the sample. This is related to the flow speed to be studied. This means that if it is high, it requires a high sampling rate to obtain reliable data, so the PRF must be high.



The WF or Wall Filter is a system to ignore low-frequency signals (and high-noise signals) from the movements of other structures, such as the vessel wall. The high filter eliminates low-speed flow underlying the base line. If no flow is detected in a vessel, it may be because the speed is very low and is being filtered. In this case, you would need to change to a low filter.

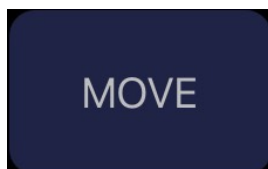

MOVE is to move the Doppler image with your finger when you have selected this button. If you press this button, it changes to SIZE.

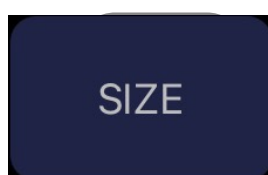

When the SIZE button is activated, you can adjust the size of the exploration area. It is advisable that the size of the "box" be as small as possible, especially as far as width is concerned, and that it be as superficial as possible, since this entails a higher number of photograms per second and the image will be formed more quickly.

3.2.3 PDI Mode

This analyses the change in amplitude and encodes it in colour scale; in other words, it represents the integrated power from the Doppler signal instead of the frequency change. Therefore, it does not provide information on flow speed, and is less dependent on the Doppler angle. It only indicates the flow presence in a structure and is more sensitive to detecting slow flow than the colour Doppler. It does not provide information on the flow direction.



PDI mode settings:

Although the analysis conducted on fluids is different than COLOUR, the settings in PDI Mode are the same as in COLOUR Mode.

See "colour mode settings" in section 3.2.2 to modify STEER, GAIN+, GAIN-, PRF, PW, MOVE and SIZE.




3.2.4 PW Mode

Also known as Spectral Doppler or Pulse Doppler, this analyses the speed of the ultra-sounded vessel flow depending on the time. The morphology of the wave shows the haemodynamic of said vessel.

PW mode settings:



If you press the STEER button, you change the tilt of the cut line, by $\pm 7^\circ$.

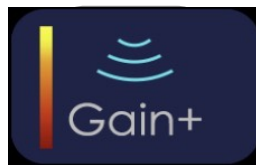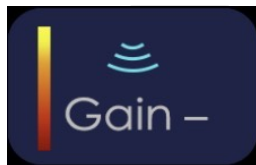
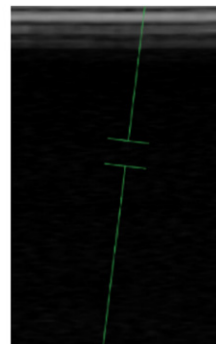


Spectral gain increases or decreases the base line signal. The most idea is that it provides for proper delineation of the wave's contour above and below.



VOLUME is the space measuring the Doppler frequency change. With Pulse Doppler, the Sample Volume is shown by a cursor (Doppler Window), whose size you can vary based on the vessel width or diameter.

It should be placed in the centre of the vessel, occupying approximately $\frac{2}{3}$ of its diameter.

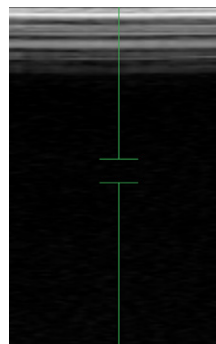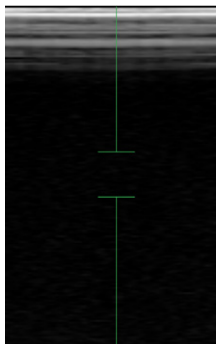




ANGLE is used to modify the measurement angle in the VOLUME window. Once the grey scale image obtained is adequate, to measure the flow speed of said structure in a reliable way, you must adjust the Doppler angle. By default, the assigned angle is 0°.

Over-correcting the angle (less than the actual angle) leads to a false low flow speed measurement, and under-correcting the angle leads to an overestimate of flow speed.




As mentioned in COLOUR Mode settings in point 3.2.2, the PRF or Pulse Repetition Frequency is a parameter set by the operator. It defines the rate or frequency with which the transducer will emit ultrasound pulses to show the sample.

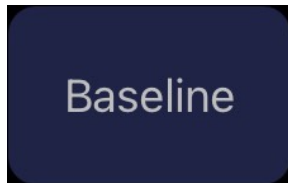

With the BASELINE function, you can modify the base line, highlighting certain aspects of the flow. When you increase flow acceleration, the PRF must be increased to prevent errors. If it cannot be increased any further and the curve is predominantly positive or negative, the scale can be adjusted by raising or lowering the base line.

Measurements in PW mode

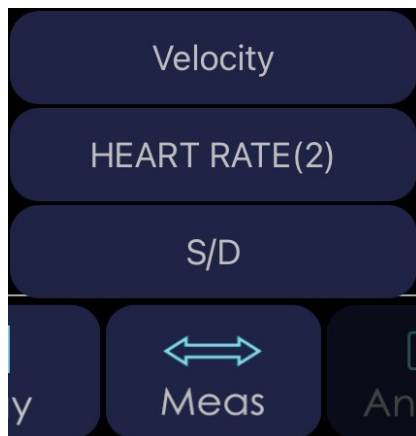


Once the examination is complete, in PW mode we can carry out the following measurements:

VELOCITY: Vertically place the cursor on the most distal part of the base line to calculate flow speed.

HEART RATE (2): By measuring with two cursors on the horizontal line and between 2 peaks, you calculate the heart rate. With the scale shown, the final result is double the actual frequency.

S/D: This is to make cardiac measurements of the Systolic Peak (PS), the End Diastolic (ED) and the relationship between the two of them (S/D).

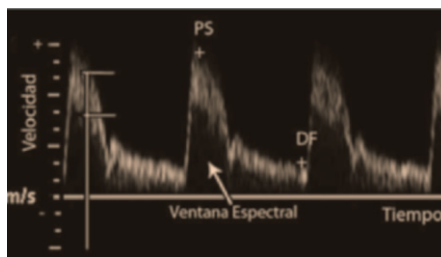

Note: All settings must be adjusted at each moment of the examination, based on the area examined, the depth of the organ to be examined, the diameter of blood vessels, flow in said vessels and the scanning angle.

4. Maintenance

4.1 Charging the probe

When the battery runs out, you must recharge it. You can do this with a wireless charger. Not included as part of the series is the loading platform. Connect its USB cable to it, and the mobile device charger.

When the charger is connected to the USB cable, the platform's LED light turns green.

We recommend using a 2A and 5V, charger and the original cable for the wireless charging base.

When the battery is being charged, the indicator levels will blink until charging is complete. At this point, the four levels will be visible, without blinking. Charging the battery can take up to 6 hours.



The battery duration is 6 hours under constant operation, and 14 hours on stand-by. As such, if you are not examining between one animal and the other, we highly recommend freezing the image to save battery if you need to work for many hours throughout the day.

The iPad lasts around 9 hours, depending on ambient light. If there is much light, the iPad's brightness will turn up to offset it, leading to higher consumption.

4.2 Tablet resistance

To prevent the tablet from becoming dirty from dust or getting moist or wet, a water-proof cover is included. It is splash- and dust-resistant.

As a particular recommendation, we suggest purchasing an IP68 cover for iPad mini, which is also fall-resistant up to 1.2 m.

4.3 Cleaning

The probe must periodically be cleaned. The part of the probe that is inserted into the rectum has IP68 water- and dust-protection, meaning it resists dust and water, but the (beige) electronic part **cannot be submerged**. You may use a damp cloth or paper to clean it.

4.4 Storage

While storing the probe, we highly recommend that the temperature be no lower than 15°C. If the probe or the mobile device's temperature is lower, you must warm them before use.

Another recommendation: regarding battery charge, if the probe and the tablet will not be used for a long period, we recommend storing them with a full charge.

5. Anexo

For the protection of the environment and in accordance with current regulations, the product and its electronic accessories should not be disposed of with other household waste at the end of its useful life. For information on the recycling of this product, please contact local authorities, the waste collection service or the manufacturer. In our firm commitment to the environment, we put at the service of the user the free collection of the equipment to discard it according to current regulations. The user will be given a certificate of this recycling process. For more information, consult your sales assistant.

6. Declaration of Conformity

NEW VETEC SL hereby declares that the type of radio equipment (probe model) is in accordance with Directive 2014/53/EU.

New VeTec S.L.	Expediente Técnico	
	Ecógrafo PROVETSCAN SR-2C	Marzo 2021
		Edición 1



DECLARACIÓN DE CONFORMIDAD UE
(de acuerdo con la norma ISO/IEC 17050-1:2004)
EU DECLARATION OF CONFORMITY
(according to ISO/IEC 17050-1:2004)

Fabricante del producto: New VeTec S.L.
Product manufacturer

New VeTec S.L.

Dirección: Calle Moisés de León nº 52 Bajo. 24006 de León.
Address

DECLARA BAJO SU RESPONSABILIDAD QUE EL PRODUCTO:
DECLARE UNDER THEIR RESPONSIBILITY THAT THE PRODUCT

Nombre: ECÓGRAFO DE ALTA RESOLUCIÓN PARA VACUNO Y EQUINO
Name: HIGH RESOLUTION ULTRASOUND SCANNER FOR DIAGNOSIS IN CATTLE AND HORSES

Marca/Modelo: New VeTec S.L. / PROVETSCAN SR-2C
Model/Reference

CUMPLE LOS REQUISITOS DEL REGLAMENTO:
CONFORMS WITH THE REQUISITES OF THE REGULATION

Directiva 2014/53/UE sobre Equipos Radioeléctricos
Directiva 2014/30/UE de Compatibilidad Electromagnética
Directiva 2011/65/UE sobre RoHS
Directive 2014/53/EU relating RED
Directive 2014/30/EU relating EMC
Directive 2011/65/EU relating RoHS

NORMATIVA APLICABLE:
APPLICABLE STANDARDS:

ETSI EN 301 489-1 V2.1.1	IEC 62321-3-1:2013	UNE-EN 55032:2015
ETSI EN 301 489-17 V3.1.1	IEC 62321-4:2013+A1:2017	UNE-EN 55024:2010
ETSI EN 300 328 V2.1.1	IEC 62321-5:2013	UNE-EN 61000-3-2:2014
UNE-EN 62479:2010	IEC 62321-6:2015	UNE-EN 61000-3-3:2013
UNE-EN 60950-1+ A11:2009 + A1:2010 + A12:2011 + A2:2013	IEC 62321-7-1:2015	EN IEC 63000:2018
	IEC 62321-7-2:2017	IEC 62321-8:2017

Marzo (March) 2021, Madrid
Fecha y lugar (Date and place)

Nombre, cargo y firma
(Name, position and signature)



New Vetec, S.L.

CIF: B42834465

Moises de León, 52 Bajo - León (España)

Tel. 633 913 164 e-mail: roberto@provetscan.com

www.newvetec.com www.provetscan.com